

HPLC REPORT


Sample: A-16
Analyst: WGJ
Lot. No.: P160106-MX630633
Column: SinoChrom ODS-BP, 4.6*250mm, 5µm
Solvent A: A: 0.1% Trifluoroacetic Acid in 100% Acetonitrile
Solvent B: B: 0.1% Trifluoroacetic Acid in 100% Water

Analyzed date: 2016-1-11

Gradient:

	A	B
0.0min	13%	87%
25.0min	38%	62%
25.1min	100%	0%
30.0min		Stop

Volume: 20µl
Wavelength: 220nm
Flow rate: 1.0ml/min

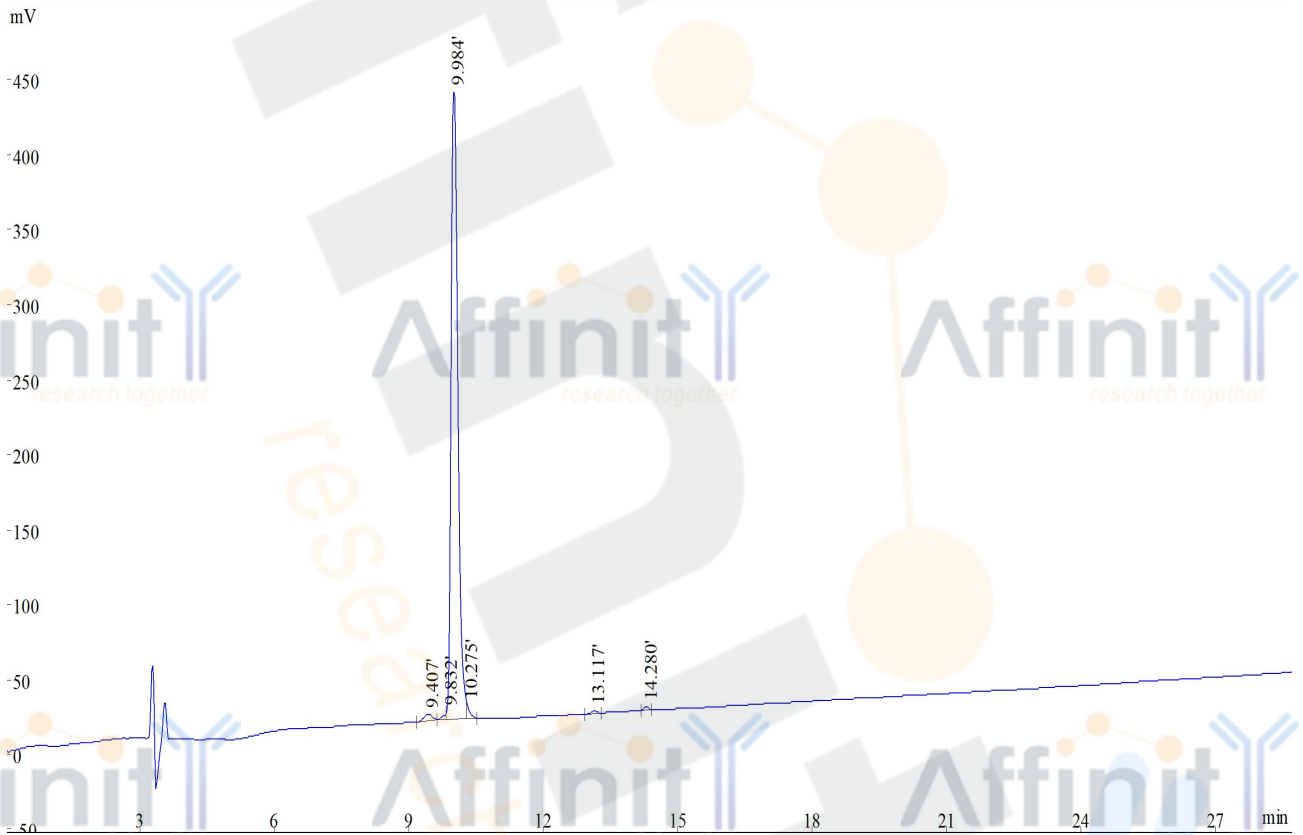

Rank	Time	Area	Conc.
1	9.407	49251	1.223
2	9.832	19863	0.4932
3	9.984	3896142	96.74
4	10.275	37168	0.9229
5	13.117	11122	0.2762
6	14.280	13734	0.341
Total			100